

ZymoPURE™


One Technology. All of your Plasmid Needs.

**25% OFF
List Price!**

Please quote CS#1977.
Offer ends
31.12.17

Innovation


Novel Binding System
up to 4x higher yield

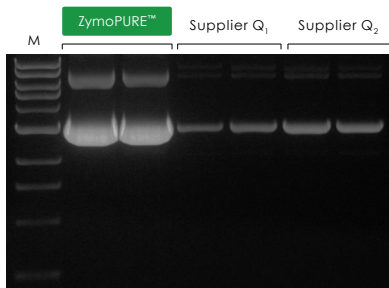

Concentrated and
Ready-to-use

Pure

EndoZero™ Column Technology
superior endotoxin removal

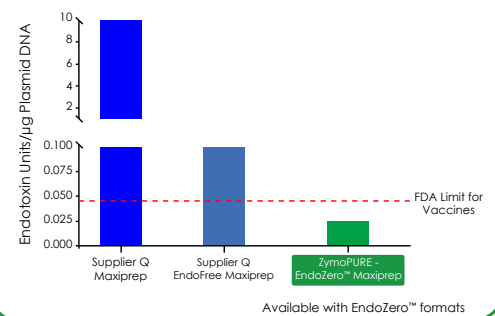

Simple

Streamlined, Pellet-free Workflow


direct lysis
no cell pelleting
required



neutralize
using Zymo's
patented color
buffer system



bind/wash
rapid loading onto
a spin-column via
vacuum or centrifuge

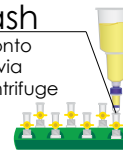

elute

transfection grade
plasmid DNA



Available with Express™ formats

Learn more at www.zymoresearch.com/zypure

	Max. Capacity	Processing Time	Culture Volume	Endotoxin Levels	Size (Cat. No.)
ZymoPURE™ Plasmid Miniprep Kit	100 µg	≤ 15 min.	5 ml	≤ 1 EU/µg	50 preps (D4209) 100 preps (D4210)
ZymoPURE™ Plasmid Midiprep Kit	300 µg	18 min.	50 ml	≤ 1 EU/µg	25 preps (D4200) 50 preps (D4201)
ZymoPURE™ Plasmid Maxiprep Kit	1200 µg	18 min.	150 ml	≤ 1 EU/µg	10 preps (D4202) 20 preps (D4203)
ZymoPURE™ Plasmid Gigaprep Kit	10 mg	40-50 min.	2.5 L	≤ 1 EU/µg	5 preps (D4204)
ZymoPURE-EndoZero™ Plasmid Maxiprep Kit	1200 µg	20 min.	150 ml	≤ 0.025 EU/µg	10 preps (D4205)
ZymoPURE-EndoZero™ Plasmid Gigaprep Kit	10 mg	50 min.	2.5 L	≤ 0.025 EU/µg	5 preps (D4207)
ZymoPURE-Express™ Plasmid Midiprep Kit	1200 µg	15 min.	25 ml	≤ 1 EU/µg	25 preps (D4213)